

T-DOC USB HID scanner configuration


To configure a T-DOC Corded Scanner (currently Datalogic QuickScan QD24XX) as a USB HID scanner the next sequence of QR-codes should be scanned one by one



Select USB Keyboard



Restore Custom Default Configuration

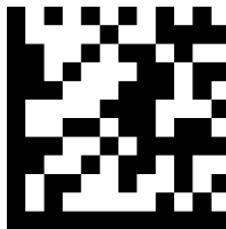

Enter Programming Mode




Enable Interleaved 2/5 (1 2 of 5)



I 2 of 5 Check Character Calculation = Check Standard (Modulo 10)



I 2 of 5 Check Character Transmission = Don't Send

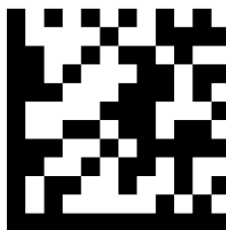

Exit Programming Mode

T-DOC standard scanner configuration

To reset/reconfigure the scanner to the standard T-DOC configuration:

1. Scan:



Select USB-COM-STD

2. Press **Program scanner** in the **Settings** dialog (activated from the **Settings** menu) in the T-DOC Scanner client
- or -
Select **Actions > Program scanners > <Scanner Name>** in the Main menu T-DOC Admin.